



Report 3

Fri Sep 08 2023 10:12:25 GMT-5


OPENSIGNALS VERSION: Public Build 2022-05-16

Author: PLUX

Devices

name				address				type				sampling rate							
98:D3:51:FD:71:01				98:D3:51:FD:71:01				bitalino_rev				1000 Hz							
channel		label		sensor		resolution		channel		label		sensor		resolution					
A1	A1	EMGBITREV		10		I1	INPUT	BTN		1		I2	INPUT	Generic Input		1			
								O1	OUTPUT	LED		1		O2	OUTPUT	BUZ		1	
								DAC	DAC	PWM		8							

Electromyographic analysis

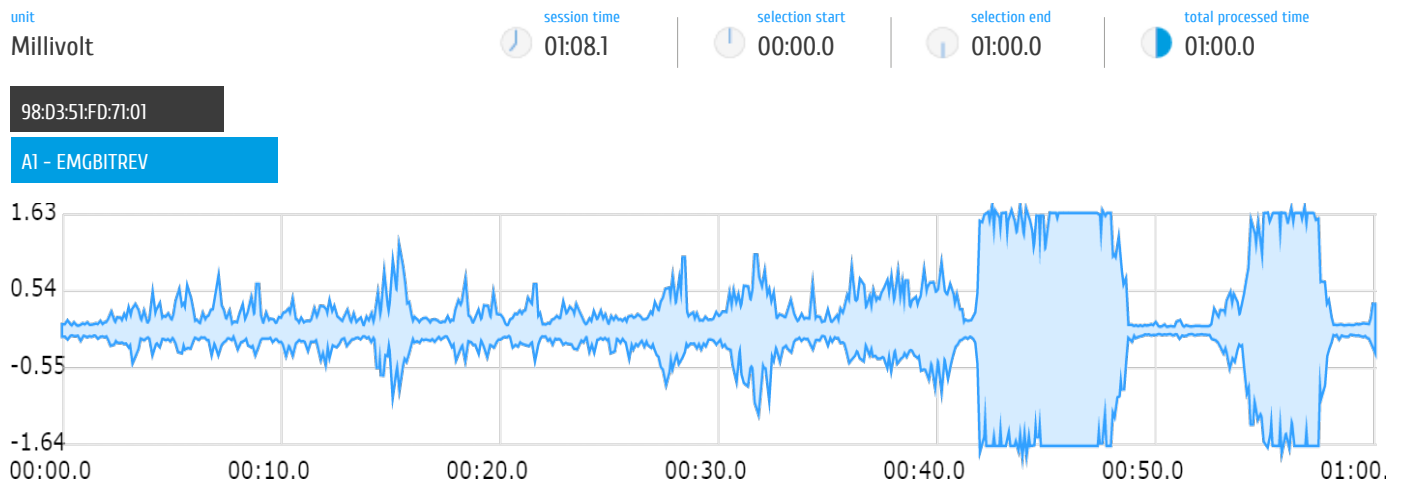

Processed channels

98:D3:51:FD:71:01

A1

FILTER
applied bandpass - [10,100] Hz cutoff - Order: 6

DETECTION OPTIONS
applied to TKE: smooth window of 0.5 seconds

Muscle Activation results

Activations

find [all muscle activations](#)

98:D3:51:FD:71:01

channel	digital	activation	start (s)	end (s)	duration (s)	latency (s)	max (mV)	min (mV)	ampl+ (mV)	ampl- (mV)	ampl+- (mV)	rms (mV)	area (mV^2)
1	NONE	1	41.889	48.309	6.420	41.889	1.63 mV	-1.63 mV	99.00	99.00	3.27 mV	0.79 mV	4.1e+3 mV^2
	NONE	2	54.203	57.546	3.343	54.203	1.63 mV	-1.63 mV	99.00	99.00	3.27 mV	0.71 mV	1.9e+3 mV^2
Diff	-	-	-	-	-	-	-	-	-	-	-	-	-
Ratio	-	-	-	-	-	-	-	-	-	-	-	-	-

Glossary

1. MVC: Maximum Voluntary Contraction.
2. TKE: Teager-Kaiser Energy operator.